



CELEBRATING 60 YEARS

Hospital



Infection Control Equipment

Utensil washer disinfectors

Stainless steel sluice units

Easy clean utensils

Stainless steel work surfaces, drainers, shelves and utensil Racks

Stainless steel deep sinks, dedicated handwash bowls, taps and traps

leaders in infection control equipment for 60 years

Stanbridge Ltd Beaumont House, Cray Valley Road, Orpington, Kent BR5 2UB

Tel: 01689 806500 email: sales@stanbridge.co.uk Web: www.stanbridge.co.uk

Stanbridge reserves the right to make design changes without notification



MAY 2015




Stanbridge specialises in sluice room and infection control equipment.

Established in the 1950s and credited with pioneering design for modern washer/disinfector machines, Stanbridge has the expertise and practical experience that can assist and reassure in every situation.

- infection control solutions for hospitals and care homes
- involvement in development of infection control procedures
- advising on sluice room requirements
- training programmes for engineers
- advice on HTM 2030 and BSEN ISO 15883-1 and -3
- participation in seminars and study days
- support and encouragement for the care sector

The comprehensive Stanbridge range includes washer/disinfectors, stainless steel products for sluice rooms, bespoke equipment and fittings, and utensils.

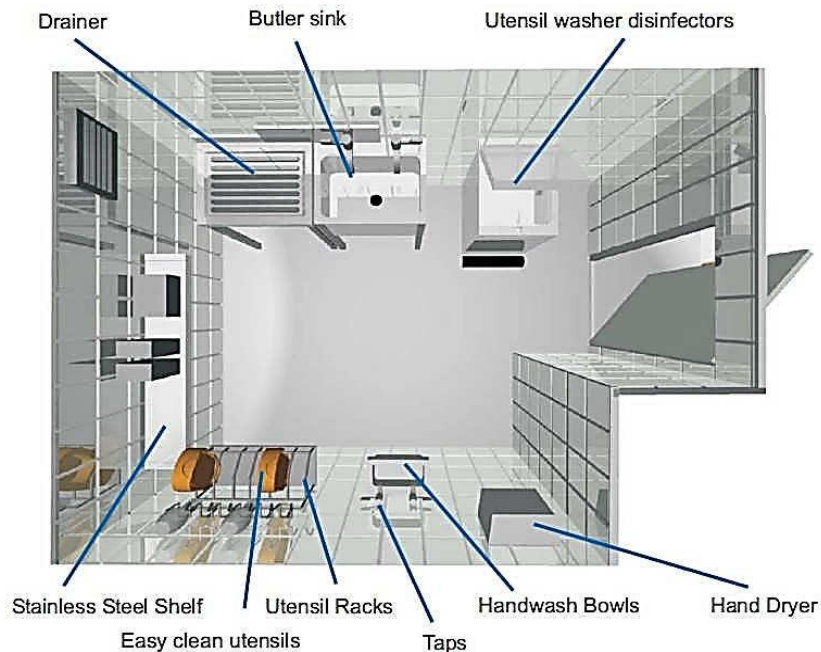

FOLLOW US ON





PLANNING A SLUICE ROOM

Preface

These notes are about controlling infection.

Guidelines and monitoring procedures are being progressively imposed to develop quality care awareness.

The control of infection covers a wide range of procedures in adult care, both domestic and Nursing where the transfer of infection between residents in group living is a natural hazard. Human waste is a significant source of this risk, but managed routines particularly in the dirty utility room will reduce any cross contamination to acceptable levels.

Introduction

We have produced this information to provide a brief guide to assist those concerned with the planning and equipping of sluice rooms.

All establishments concerned with care, regardless of size, have to provide facilities to ensure the hygienic disposal of human waste.

Planning

Ideally a sluice room should be planned at the outset just like the kitchen and the laundry. The number of sluice rooms will depend on the number of floors and the layout of the building. The functions and the services required will have to be taken into consideration, such as:

Hot and cold water inlets

Power supplies and voltages

Soil outlets and vents

Floor and wall finishes

Adequate ventilation

Remember that a sluice room can be used 24 hours a day, so a site should be chosen where noise from water or clattering hardware is not a nuisance.

Wherever possible the room should be divided into 'clean' and 'dirty' areas. Each room needs a sink, preferably with a drainer. Taps are a source of contact by hands and automatic taps or those with elbow levers are preferred. Push button or foot operated soap dispensers are recommended.

Dedicated hand washing facilities are needed.

Equipment

The term 'Sluice Room' is used loosely to define a utensil washer/disinfector with a back up of a WC and a designated sink for washing.

The automatic washer/disinfector for commode bucket, bedpan and urinals is, however, a far more sophisticated piece of hardware than has long been the norm in NHS hospitals. The new generation of machines can carry out the functions of both slop hopper and washer/disinfector.

The term 'disinfection' is used to describe the 'bug kill' cycle. Moist heat is used from an incorporated steam generator to heat each utensil. It is NOT liquid from a bottle.

Our machines have been engineered over years of testing, establishing effective cleaning and disinfecting of the specific utensils used in human waste disposal.

Designed for a function, normally they will accommodate bedpans, commode buckets, urine bottles or various bowls/utensils. The machines open with a pneumatic foot pedal to give the operator full hand control of the receptacle. They have cold and hot washes, the cold wash removes the soil, the hot wash removes extra stubborn waste, then items are disinfected. The minimum requirement is 80 deg. C for 60 seconds. The machines usually attain 90 deg. C+ for 70 seconds, followed by a cold rinse to cool the receptacle. Interlocks prevent opening during the cycle. Utensils removed from the machine should be placed on racks adjacent to the machine - this is the clean area.

If used soiled items have to be left prior to washing they should have a separate 'dwelling place' - not in the clean area. A trolley or stainless steel worktop should be provided.



The SLIMLINE HS WASHER/DISINFECTOR



The new generation of Stanbridge Washer/disinfectors has been designed to improve infection control on today's busy and challenging Hospital wards.

One Machine to –

- ✓ Dispose of utensil contents
- ✓ Thoroughly clean and remove soil
- ✓ Thermally disinfect by steam to remove pathogens and micro organisms ensuring safe reuse.

Features Include –

- A unique vertical rotating column jet system for extensive and safe cleaning.
- Newly designed load carrier to accommodate today's range of utensils.
- Digital programmed display to visually prove disinfection.
- Hands free lid/door opening operation.
- Descalant and detergent system fitted as standard.
- Cycle counter to assist the organization of servicing and validation to ensure stringent requirement intervals are met.

Additional Safety Features –

- Automatic cut off for blockages.
- Automatic lid/door lock during the cycle operation.

For ultimate infection control the rigid support carrier for cardboard liner bed pans must be treated as a reusable utensil.

For your peace of mind our range of washer/disinfectors provide the solution to all user needs for re-usable items which include pulp support carriers ensuring clean and disinfected utensils every time.



Size
Height
Width
Depth

Installation
 1500mm
 535mm
 490mm

Electrical

230V 50Hz 16A single phase
 Fuse - mains 16A - control 5A

Water Feed
Supply
Connection

Cold Feed x 2 Via Solenoid Valve
 1/2" BSP male parallel

Washing
cycles

1. Bottle Wash
 Wash - Disinfection at 80°C for 70 seconds -
 Rinse to cool
2. Pan Wash
 Pre flush - Wash -
 Disinfection at 80°C for 70 seconds - Rinse to cool

Overflow

2 x 45mm x 30mm Rectangle Slot
 End must be visible BS EN 15883-1
 (UK requirements).

Waste

82mm (3") 'P' or 'S' trap



Stanbridge Limited, Beaumont House, Cray Valley Road,
 St. Mary Cray, Orpington, Kent BR5 2UB UK
 Tel +44 (0) 1689 806500 www.stanbridge.co.uk





The CS1 HS WASHER / DISINFECTOR



The Compact Superior 1 HS washes and disinfects a variety of utensils one at a time and is specially designed for use where space is limited.



The CS1 provides an effective infection control barrier which protects residents and staff from the risk of cross infection and management from potential liabilities. All Stanbridge washer/disinfectors include the DEOSCALE water treatment system to control limescale build-up and reduce maintenance

The new generation of Stanbridge Washer/disinfectors has been designed to improve infection control on today's busy and challenging Hospital wards.

Features Include –

- A unique vertical rotating column jet system for extensive and safe cleaning.
- Newly designed load carrier to accommodate today's range of utensils.
- Digital programmed display to visually prove disinfection.
- Hands free lid/door opening operation.
- Descalant and detergent system fitted as standard.
- Cycle counter to assist the organization of servicing and validation to ensure stringent requirement intervals are met.

Efficient tested results achieved due to the fixed and rotating nozzle/candle which provides optimum water distribution inside and outside the utensil.

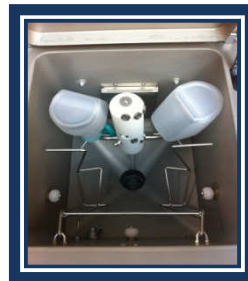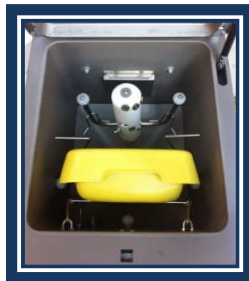

<p>Size</p> <p>Height</p> <p>Width</p> <p>Depth</p>	<p>Installation</p> <p>915mm</p> <p>435mm</p> <p>650mm</p>	<p>Electrical</p> <p>230V 50Hz 16A single phase</p> <p>Fuse - mains 16A - control 5A</p>	
<p>Water Feed Supply Connection</p>	<p>Cold Feed x 2 Via Solenoid Valve</p> <p>1/2" BSP male parallel</p>	<p>Washing cycles</p>	<p>1. Bottle Wash</p> <p>Wash - Disinfection at 80°C for 70 seconds -</p> <p>Rinse to cool</p> <p>2. Pan Wash</p> <p>Pre flush - Wash -</p> <p>Disinfection at 80°C for 70 seconds - Rinse to cool</p>
<p>Overflow</p>	<p>45mm x 30mm Rectangle Slot</p> <p>End must be visible BS EN 15883-1 (UK requirements).</p>		
<p>Waste</p>	<p>82mm (3") 'P' or 'S' trap</p>		



Stanbridge Limited, Beaumont House, Cray Valley Road,
St. Mary Cray, Orpington, Kent BR5 2UB UK
Tel +44 (0) 1689 806500 www.stanbridge.co.uk





The CS3 HD HS WASHER / DISINFECTOR



The Compact Superior 3 HD HS **multi-utensil** machine.

The CS3 HD HS has the capacity to efficiently wash and disinfect 3 utensils in one operation, cutting water and heating costs, staff waiting time and utensil queueing. Stanbridge continues to increase its reputation through its investment in product development.

This range of Washer/disinfectors are consistently tested and efficient washing performance.

Efficient tested results achieved due to the fixed and rotating nozzle/candle which provides optimum water distribution inside and outside the bedpan/utensil.


Steam is distributed through the candle to ensure the complete disinfection of the utensil during each cycle.

To avoid dirt and bacteria growth, the internal part of the chamber is manufactured with smooth edges

Additional Safety Features –

- Automatic cut off for blockages.
- Automatic lid/door lock during the cycle operation.

For ultimate infection control the rigid support carrier for cardboard liner bed pans must be treated as a reusable utensil.

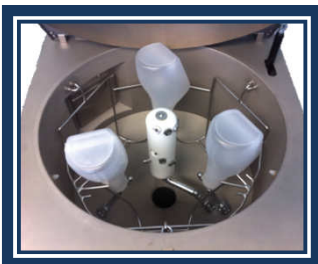

Size
Height
Width
Depth

Installation
1100mm
600mm
712mm

Electrical

230V 50Hz 16A single phase
Fuse - mains 16A - control 5A

Water Feed
Supply
Connection

Cold Feed x 2 Via Solenoid Valve
1/2" BSP male parallel

Washing
cycles

1. Bottle Wash
Wash - Disinfection at 80°C for 70 seconds - Rinse to cool

Overflow

45mm x 30mm Rectangle Slot
End must be visible BS EN 15883-1
(UK requirements).

2. Pan Wash

Pre flush - Wash - Disinfection
at 80°C for 70 seconds -
Rinse to cool

Waste

82mm (3") 'P' or 'S' trap



Stanbridge Limited, Beaumont House, Cray Valley Road,
St. Mary Cray, Orpington, Kent BR5 2UB UK
Tel +44 (0) 1689 806500 www.stanbridge.co.uk





PULP SUPPORT WASHER/DISINFECTOR



For ultimate infection control the rigid support carrier for cardboard liner bed pans must be treated as a reusable utensil.

The Stanbridge range of washer/disinfectors are perfect to ensure reusable carriers are washed and disinfected to the standards detailed in HTM 20/30 and BS EN ISO 15883.

For your peace of mind the Stanbridge range of washer/disinfectors provide the solution to all user needs for reusable products, bedpans, commode pots, bowls, bottles, jugs and support carriers ensuring clean and disinfected utensils every time.





Water treatment by DeoScale

DeoScale is a water additive to aid cleaning and control the build-up of hard scale deposits in automatic machines that wash and disinfect reusable bedpans and commode pots.

Area Scale

In some areas the build-up of limescale is high enough to reduce the effectiveness of the washer/disinfector by causing surface deterioration of utensils. Both performance and appearance suffer.

Why water treatment?

Water for washing containers for human waste should be conditioned in the same way as dishwashers for household plates and cups to ensure the fatty residue is removed by reduction in the water surface tension.

The centre of excellence for clinical safety of utensils is hosted by the Department of Health and contained in the recommendations of HTM20/30 for best practice.

Guidelines for the use of washer/disinfectors for handling human waste containers include necessary wash treatments to achieve good outcome.

DeoScale fluid

DeoScale fluid has two components -a non-foaming detergent and a descaler. The addition of a small measured dose of DeoScale fluid to the wash cycle reduces hard water scale build-up. DeoScale is a low-hazard liquid detergent/descaler designed for auto dosing. It is efficient in removing organic and inorganic soils. Low dose in soft water areas reduces water browning.

Background

There are many 'fixed gadget' devices offered to control or treat limescale in water supplies. Confidence in outcome is inconsistent usually because of a dependency on flow rate and limescale levels.

In consequence DeoScale has been developed to control scale build-up in water systems as a new technology for bedpan or commode washer to provide a predictable outcome for varying degrees of district hardness as well as water flow rates. The dose can be set to match local conditions.

Benefits

The system is based upon the accurate metered introduction of the specially formulated fluid into one break tank during each cycle. The chemical composition of this fluid is such that the component of temporary hardness in the water is disorientated and prevented from forming a hard scale layer. Any precipitate passes to drain in each water flow cycle.

New equipment remains pristine. Scaled chambers in old equipment return to 'as new' after a once only hand descale. Existing scale on utensils over a period slowly re-combines and surfaces to 'bright'. Pipe work is cleaned. Organic grubbiness is removed.

Breakdowns and engineers' time has fallen on trial sites as the internal system scale reduces by slowly returning to a soft flushable character.

Cost effectiveness can be demonstrated. Utensil management is excellent and in consequence, so is care staff reaction. Soil removal and control of infection is enhanced.

Based upon the water provider's analysis of the water supplied to the washer/disinfector, the dose required can be calculated or put to work on a 'trial and error' basis.

The count of machine cycles per day for each dirty utility room will identify a budget cost.

Hazard is low - liquid is classified as a mild Irritant

Equipment

The basic dosing facility can be fitted in the factory for new installations or site added and consists of:

1. Metering pump
2. Fluid carrier for a 5 litre container suitable for wall mounting

DeoScale fluid is packed in 5 litre screw top plastic containers.

Typical sales pack is 4x5 litre boxed suitable for groupage transport.

Contract quantities per year are negotiable. Please contact our Spares Department for prices.

Consumption

The DeoScale consumption for a typical care home or ward with 30 beds would be around 1x5 litres per month for a single machine. In a high scale district this may rise to 2x5 litre per month.



Stanbridge Limited, Beaumont House, Cray Valley Road,
St. Mary Cray, Orpington, Kent BR5 2UB UK
Tel +44 (0) 1689 806500 www.stanbridge.co.uk





POTS - PANS & BOTTLES



Commode Insert pan

The insert commode pan improves the commode chair use, as well as sample gathering in the toilet suite.



Slipper pan

The slipper pan has no end pocket to trap soil and eliminates cleaning problems. Three hand grips allow easy positioning. An indent on the plateau reduces point load on the spine. Smooth movement over the mattress is assured.



In-Bed pan

The in-bed pan gives confidence, stability and secure feeling for sitting patients. The buttress edge gives extra support on the mattress, especially backwards. A raised edge helps location and comfort.



Urine bottle


Robust opaque polythene urine bottle, graduated in ml to 1 litre. Complete with snap lid for transit safety.



PAN BOTTLE & UTENSIL RACK

The Amethyst Economic (non-heated) quality stainless steel racks are designed to store bedpans, commode pots and urine bottles efficiently and safely.

All racks are manufactured to an exceptionally high standard in electro polished grade 304 stainless steel with a removable stainless drip tray for easy cleaning and disinfection.

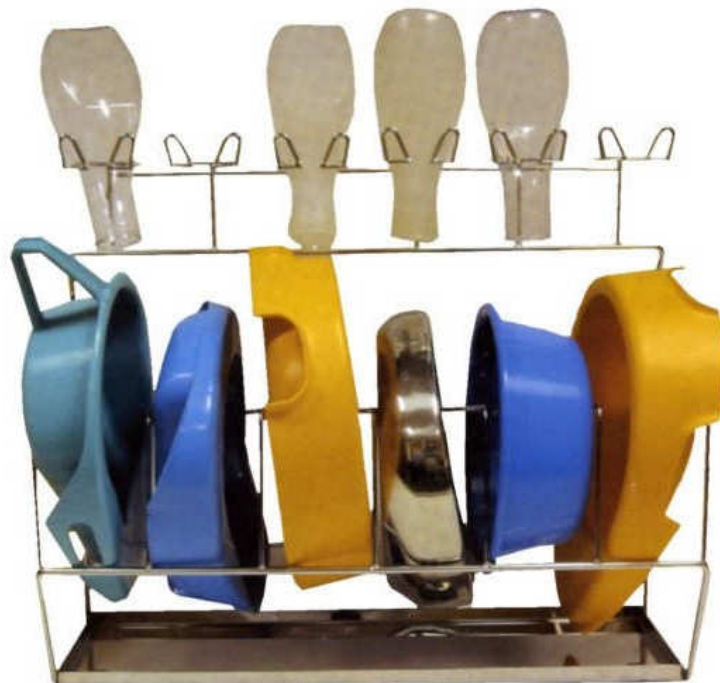



HAND WASH BASINS



WHB320R

Hand basin constructed of grade 304 stainless steel. Finish bright/satin. Complete with skirt/apron, bolt set, waste fitting /NO overflow and plug. Taps and traps extra.

Bowl size 320dia x 105mm depth

Overall size: W400 x D350 x H105mm



STHB02

Corner mounted hand basin constructed of 0.9mm thick 304 grade satin finish stainless steel for corner fixing and complete with waste. Taps and traps extra.

Bowl size: 212 dia x 90mm deep.

Overall size: W305 x D305 x H90mm



STAN5

Wall mounted hand basin manufactured from 1.2mm thick 304 grade bright/satin finish stainless steel. Tap holes 28mm dia.

Supplied with skirt, wall fixing.

Taps and traps extra.

Bowl size: W300 x D150 x H100mm

Overall size: W350 x D260 x H185mm



JANITORIAL EQUIPMENT & DEEP CLEAN SINKS

Stainless steel deep sinks, dedicated handwash basins, taps and traps



Combined Janitorial & Hand Wash Unit

Manufactured from 304 grade stainless steel and supplied with adjustable legs, hinged bucket grid, bowl with monobloc swivel mixer tap and combined waste and overflow.

Dimensions: W500 x D750 x H880mm



Free Standing & Wall Mounted Butler Sinks

Manufactured from 1.2mm thick 304 grade satin finish stainless steel with 38mm waste, 2 tapholes, NO overflow or plug, 100mm rear upstand and anti-drip edges all round. Mounted on 30mm square tubular stainless steel uprights with adjustable feet but no undershelf. Can be used in conjunction with the Stanbridge range of stainless steel drainers or independently.

Sink 1

Bowl capacity 48 litres W400 x D400 x H300
Overall dimensions: W500 x D600 x H850



Sink 2

Bowl capacity 60 litres W500 x D400 x H300
Overall dimensions: W600 x D600 x H850

Sink 3

Bowl capacity 72 litres W600 x D400 x H300
Overall dimensions: W700 x D600 x H850



SHELVING DRAINERS & WORK SURFACES



Island work surface - Free standing

Manufactured from 304 grade stainless steel with 60mm turn down all round and a standard working height of 850mm. Supplied with stand adjustable feet, undershelf and earth tag.

Dimensions: Up to 3500mm long, width 650mm



Sink unit/Work surface

Manufactured from 304 grade stainless steel with 60mm turn down all round and a standard working height of 850mm. Supplied with 'single bowl' and rear upstand.

Dimensions: Up to 3000mm long, width 650mm
Bowl size: 610 x 457 x 300mm deep.



Wall mounted shelf/drainer

Manufactured from 1.2mm thick 304 grade satin polished stainless steel. Supplied with 2 wall mounting brackets, side aprons, 60mm rear upstand, 2 x 30mm square front support legs with adjustable feet. Butler sink to either lip over left or right.

Drainer 600-L600mm x W600mm x H960mm
Drainer 900-L900mm x W600mm x H960mm
Drainer 1200-L1200mm x W600mm x H960mm



Heavy duty wall mounted shelves

With a standard width of 300mm with integral wall fixing brackets.

Shelving 600
Shelving 1200
Shelving 1500
Shelving 1800
Shelving 2000



TAPS & BOTTLE TRAPS



STAN EC 04

Mono Bloc mixer tap sink mounted.
Chrome plated quarter turn compression valve .
3" lever length (6" lever also available)



STAN ST 04

'H' Section mixer taps sink mounted.
Chrome plated quarter turn ceramic. 6" lever.



STAN EC 03

High neck tap sink mounted.
Chrome plated quarter turn ceramic valve
3" lever length. (6" lever also available)



STAN ST 02

Bib taps for wall mounting. Chrome plated
quarter turn ceramic valve 3" lever length
(6" lever also available)



BOTTLE TRAPS

Chrome plated Bottle Traps



STAN EC 01

Hand Basin taps. Chrome plated
quarter turn ceramic valve 3" lever length
(6" lever also available)



TOPAZ PLASTER SINK



TOPAZ
(Plaster Sink and Work Surface)

Stainless steel left or right handed plaster sink and work surface.

250mm diameter stainless steel plaster sump and removable mesh filter basket with lifting handle and rebated cover plate.

38mm waste outlet and chrome plated 38mm waste fitting.

Upstand to rear and anti-drip edges to front and sides.

Chrome plated bottle trap.

Elbow lever action chrome plated mixer taps with single swivel nozzle for surface (tap holes on 200mm centres) or wall mounting.

Two sets of bearers and stainless steel front support legs with adjustable feet or two cantilever wall brackets. Earthing terminal.



MANUAL SLUICES - SLOP HOPPERS



GARNET SLUICE (Cleaner, Sluice/Sink)

Stainless steel cleaner sluice/sink unit. Stainless steel hinged bucket grid.

Stainless steel wall mounted 9 litre cistern with cover, chrome plated lever action flush handle, side entry supply and overflow and standard cistern kit.

Stainless steel flush pipe.
Upstand to rear and anti-drip edges to front and sides.

Wall mounted wrist lever action hot and cold chrome plated bib taps extended for concealed or exposed supply.

Two cantilever wall brackets.

Hopper 530 x 530mm.



SAPPHIRE SLUICE (Sluice Unit)

Stainless Steel Sluice Unit.

Stainless steel wall mounted 9 litre cistern with cover, chrome plated lever action flush handle, side entry supply and overflow and standard cistern kit.

Stainless steel flush pipe.
Upstand to rear and anti-drip edges to front and sides.

Surface mounted wrist lever action chrome plated monobloc dual flow mixer tap with swivel spout to sluice.

Two cantilever wall brackets.

Hopper 530 x 530mm.



COMBIND MANUAL SLUICES



AQUAMARINE
(Sluice and Work Surface)

Stainless steel left or right handed sluice and work surface.

Stainless steel wall mounted 9 litre cistern with cover, chrome plated lever action flush handle, side entry supply and overflow and standard cistern kit.

Stainless steel flush pipe.

Upstand to rear and anti-drip edges to front and sides.

Surface mounted wrist lever action chrome plated monobloc dual flow mixer tap with swivel spout to sluice & Bottle trap.

Two sets of bearers and stainless steel front support legs with adjustable feet or two cantilever wall brackets.

Hopper and work surface 1070x530mm.



COMBIND MANUAL SLUICES



OPAL SLUICE

(Sluice Sink and Work Surface)

Stainless steel left or right handed sluice and sink unit with work surface.

Stainless steel wall mounted 9 litre cistern with cover, chrome plated lever action flush handle, side entry supply and overflow and standard cistern kit.
Stainless steel flush pipe.

Upstand to rear and anti-drip edges to front and sides.

Chrome plated 38mm deep seated waste fitting to sink and plastic overflow pipe.

Wrist lever action hot and cold chrome plated monobloc dual flow mixer tap with swivel spout to sluice and Bottle trap.

Two sets of bearers and two stainless steel front support legs with adjustable feet or two cantilever wall brackets.

Hopper, sink and work surface 1600x530mm.



CUPBOARDS & SINK UNITS

Cupboard and Sink Unit



Constructed of AISI 304 18/8-10 stainless steel.

Available in 12 models with 2 different basin sizes to cover the vast majority of needs that may arise.

Equipped with sloping groove drainers, so as to guide the draining water into the basin.

Equipped with stainless steel legs which are adjustable in height.



Wall Mounted Cupboard

Constructed of AISI 304 18/8-10 stainless steel.

Available in 4 models of varying sizes

Equipped with one shelf which is adjustable in height.



CELEBRATING 60 YEARS

Stanbridge Limited,
Beaumont House, Cray Valley Road,
St. Mary Cray, Orpington, Kent BR5 2UB UK

+44 (0) 1689 806500 www.stanbridge.co.uk

leaders in infection control equipment for 60 years

