



The Swing Reclining Bath System



Stanbridge Ltd

Leaders in Infection Control Equipment for over 50 Years

Beaumont House, Cray Valley Road, Orpington, Kent BR5 2UB

Tel: 01689 806500 Fax: 01689 806501 email: sales@stanbridge.co.uk www.stanbridge.co.uk

Stanbridge reserves the right to make design changes without notification

B0400

THE SWING BATH

The 'Swing' reclining bath is for truly assisted bathing.

The new 'Swing' reclining bath with side access allows the new priority of the personal care package to be developed for individual service users who still have some mobility and do not need hoisting. Wheelchair as well as walking access raises self motivation and personal satisfaction. Bathing time and carer attention can be significantly reduced with direct savings - client fatigue and hence stress is lower.

Easy, accessible and safe.

The 'Swing' is a freestanding reclining bath with carer access from three sides. The bath foot well is pre-filled with thermostatically mixed water to achieve immediate user comfort when entering the bath.

Access from the opened large door allows both hoisted and non-hoisted entrance thus ensuring the safe and effective transfer whilst maintaining the bather's well being.

When reclined the water gently moves over the body and allows manual top up.

A pivoted 'shepherd's crook' hand spray mounting ensures the spray head cannot enter the water, ensuring that the type 'A' air gap is maintained, satisfying water authorities.

With optional air spa, disinfection, height variation and auto fill, the 'Swing' offers a unique bathing satisfaction to all categories of user.

Standard equipment:

- Electric powered reclining bathing system on freestanding chassis
- Thermostatic mixer scald protected, shower with adjustable spray and flexible hose
- Head cushion
- Two handgrips

Options:

- Air spa with 10 nozzles and pneumatic control box
- Height adjustment
- Full integrated disinfection system for bath and air spa
- Automatic bath fill

We have a policy of continuous product improvement and changes may occur without notification.



Stanbridge Ltd

Leaders in Infection Control Equipment for over 50 Years

Beaumont House, Cray Valley Road, Orpington, Kent BR5 2UB


Tel: 01689 806500 Fax: 01689 806501 email: sales@stanbridge.co.uk www.stanbridge.co.uk

Stanbridge reserves the right to make design changes without notification

B0401

THE SWING BATH

Space requirements

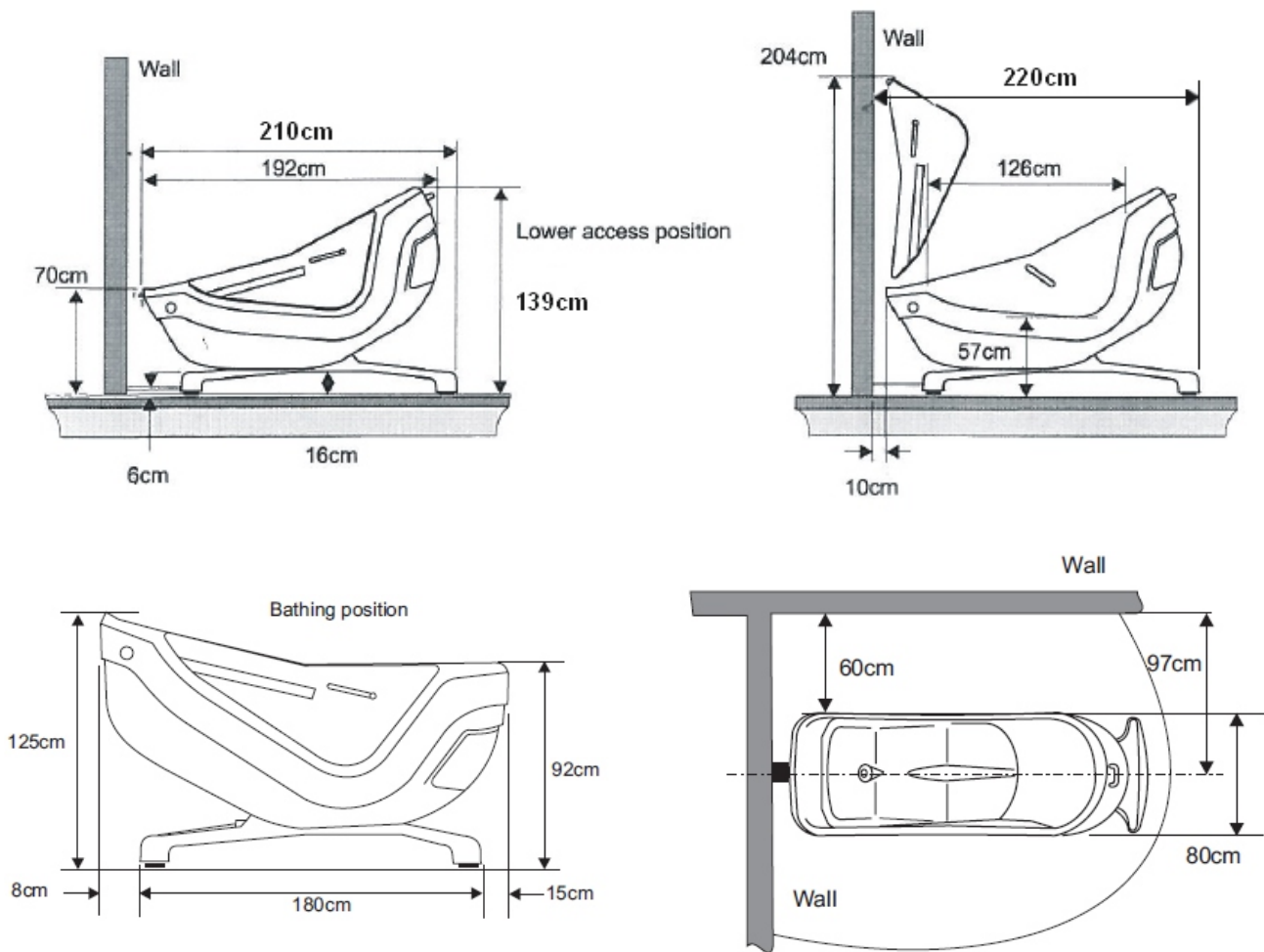

Plumbing

The Swing reclining bathtub is supplied with two flexible hoses. They must be connected to the cold and hot water supply pipes each of 3/4". Both hot water and cold water pipes should be equipped with a stop valve and a non-return valve (not supplied). The bathtub is further supplied with a flexible drainage hose with a diameter of 40mm. It should be connected to a PVC drainage pipe.

Electricity

The Swing bath system is supplied with a double insulated power cord. It should be connected without interruption, directly to the wall connection box (IP54)

Power supply: 220V/240V 50/60Hz (5A), protected by a 30mA circuit breaker.



Stanbridge Ltd

Leaders in Infection Control Equipment for over 50 Years

Beaumont House, Cray Valley Road, Orpington, Kent BR5 2UB

Tel: 01689 806500 Fax: 01689 806501 email: sales@stanbridge.co.uk www.stanbridge.co.uk

Stanbridge reserves the right to make design changes without notification

B0402