


Microtron

Made for microtitre plate applications



- ▶ Penta Drive
- ▶ Vibrationless
- ▶ High Speed
- ▶ Stackable
- ▶ Well configured
- ▶ High capacity

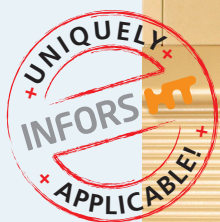

www.infors-ht.com

May be subject to technical amendments

INFORS **HT**

High capacity and high speed

▶ Penta-drive

The shaking tray is firmly connected to its support with 5 points of attachment for robust, reliable operation and minimal wear.

▶ Vibrationless

A fully optimised drive system allows high shaker speeds while ensuring unwanted noise and vibrations are eliminated.

▶ High speed

The combination of high speed and short throw make the Microtron well suited to many applications involving microtitre plates.

▶ Stackable



The Microtron can be stacked as one deck of a general Multitron system or three chambers used as a dedicated 3-stack system.

▶ Well configured

The Microtron can be optionally equipped with humidity control and/or control of CO₂.

▶ High capacity

Special platforms are available which allow 72 or even 96 x standard microtitre plates to be stacked in columns - or 48 deep well plates.

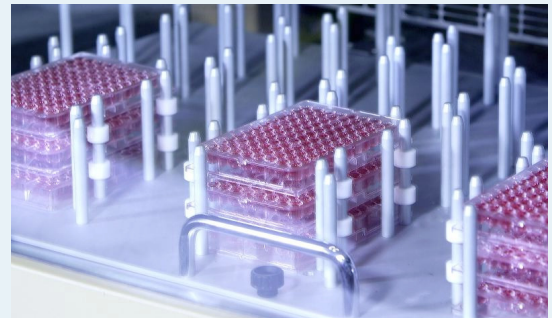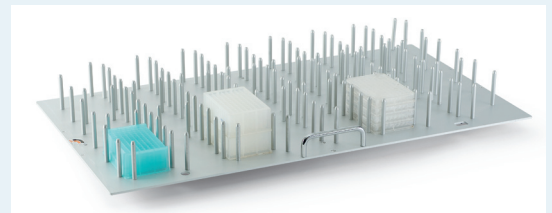

Technical Data

Shaker Speed	50–999 RPM
Throw	∅ 3 mm
Accuracy	± 1% of maximum speed
Temperature range	5 °C above RT to 60 °C
– with cooling	Approx. 15 °C under RT
Temperature resolution	0.1 °C, at the
Precision (Pt100)	± 0.2 °C
Weight without table	approx. 90 kg
Dimensions W x D x H	107 x 86 x 55 cm
Power supply	115/230V±10%, 50/60Hz

Initial capacity

Up to 72 Microtiter
– 48 deep well plates

Options

- Regulation of humidity
- Regulation of CO₂ gas
- Validation

Infors AG
Headoffice, Switzerland

Rittergasse 27
CH-4103 Bottmingen
T +41 (0)61 425 77 00
F +41 (0)61 425 77 01
headoffice@infors-ht.com

Infors GmbH

Dachauer Str. 6
D-85254 Einsbach
T +49 (0)8135/83 33
F +49 (0)8135/83 20
infors.de@infors-ht.com

Infors UK Ltd

The Courtyard Business
Centre, Dovers Farm
Lonesome Lane, Reigate
Surrey, UK-RH2 7QT
T +44 (0)1737 22 31 00
F +44 (0)1737 24 72 13
infors.uk@infors-ht.com

Infors Sarl

6, rue Marcel Paul
FR-91742 Massy Cedex
T +33 (0)169 30 95 04
F +33 (0)169 30 95 05
infors.fr@infors-ht.com

Infors Benelux BV

Logistiekweg 20-02
7007 CJ Doetinchem
P.O. Box 613
NL-7000 AP Doetinchem
T + 31(0)314-36 44 50
F + 31(0)314-37 80 76
infors.bnl@infors-ht.com

Infors Canada

8350 rue Bombardier
Anjou, Quebec
Canada H1J 1A6
T + 1 514 352 50 95
F + 1 514 352 56 10
infors.ca@infors-ht.com

**Infors Bio-Technology
(Beijing) Co., Ltd.**

Room 505C, Building 106
Lize Zhongyuan
Wangjing New Industrial Zone
Chaoyang District, Beijing
100102 P.R. of China
T +86 10 51652068
F +86 10 64390585
info@infors-ht.com.cn