

Techfors-S

The solid basis for your Pilot Fermentation



- ▶ All around access
- ▶ Compact
- ▶ Easy navigation
- ▶ Simple yet Flexible
- ▶ Low volume Specialist



www.infors-ht.com

May be subject to technical amendments

INFORS HT

High-Tech for an entry-level price

▶ All around access

Techfors-S has the inlet and exit filters mounted behind the top plate surface for easy access to all parts. Being mobile, Techfors-S can even be swivelled for rapid access to the open-frame components at the rear.

▶ Compact

The single column design of Techfors-S makes the unit compact and capable of fitting into almost any laboratory. A steam generator can be integrated into the same footprint, if required.

▶ Easy navigation

The Touchscreen ontroller uses tabs and a colour touch screen display for easy, intuitive navigation. Both new users and advanced operators are catered for.


*The X-DDC operating panel with navigation knob is still available as an option on Techfors and Techfors-S – see separate brochure for details.

▶ Simple yet Flexible

Techfors-S was designed to bring in-situ sterilizable systems to people new to this class of fermenter and it achieves this while allowing for expansion, re-configuration and open-frame gas mixing.

▶ Low volume Specialist

The jacket of each Techfors-S vessel extends below the sides of the vessel, providing good temperature control even when low volumes are used - vital for modern fed-batch applications.

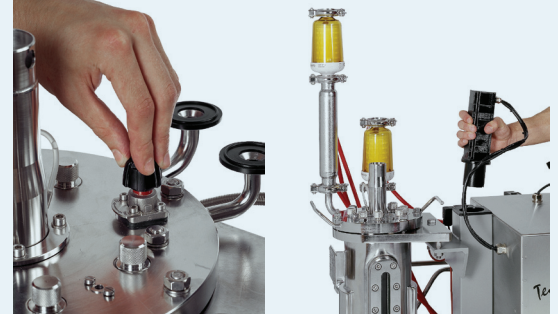

Specification

Bacterial and Cell Version:

• Vessel size	7,5 L ; 15 L or 30 L TV
• Sterilisation	In Situ
• Temperature	Ambiant/coolant +5 - 70°C
• Vessel types	316L; Double Jacket
• Speed	
Bacterial Version	50–1200 RPM
Cell Version	20–300 RPM
• Vessel Seal	
Bacterial Version	Mechanical seal
Cell Version	Magnetic coupling
• Gassing options	
Bacterial Version	Rotameter (options: mass flow valve and oxygene supplementation)
Cell Version	Rotameter and gas mix (options: mass flow valve)
• Number of pump	3 (option: +1 pump)
• Parameters	°C, RPM, pH, pO ₂ , AF, gas flow, level, pressure, weight, feed, analog pump

Fly001_en_1208

Infors AG
Headoffice, Switzerland

Rittergasse 27
CH-4103 Bottmingen
T +41 (0)61 425 77 00
F +41 (0)61 425 77 01
headoffice@infors-ht.com

Infors GmbH

Dachauer Str. 6
D-85254 Einsbach
T +49 (0)8135/83 33
F +49 (0)8135/83 20
infors.de@infors-ht.com

Infors UK Ltd

The Courtyard Business
Centre, Dovers Farm
Lonesome Lane, Reigate
Surrey, UK-RH2 7QT
T +44 (0)1737 22 31 00
F +44 (0)1737 24 72 13
infors.uk@infors-ht.com

Infors Sarl

6, rue Marcel Paul
FR-91742 Massy Cedex
T +33 (0)169 30 95 04
F +33 (0)169 30 95 05
infors.fr@infors-ht.com

Infors Benelux BV

Logistiekweg 20-02
7007 CJ Doetinchem
P.O. Box 613
NL-7000 AP Doetinchem
T + 31(0)314-36 44 50
F + 31(0)314-37 80 76
infors.bnl@infors-ht.com

Infors Canada

8350 rue Bombardier
Anjou, Quebec
Canada H1J 1A6
T + 1 514 352 50 95
F + 1 514 352 56 10
infors.ca@infors-ht.com

Infors Bio-Technology
(Beijing) Co., Ltd.

Room 505C, Building 106
Lize Zhongyuan
Wangjing New Industrial Zone
Chaoyang District, Beijing
100102 P.R. of China
T +86 10 51652068
F +86 10 64390585
info@infors-ht.com.cn