

# BTL-6000 LYMPHASTIM

BTL PRESSO-THERAPY SYSTEM FOR MEDICAL AND AESTHETIC PRACTICES



## TECHNOLOGY

The human lymphatic system is responsible for removing interstitial fluid from tissues. It absorbs and transports fatty acids and fats throughout the circulatory system and transports antigen-presenting cells to the lymph nodes where the immune response is stimulated. Diseases and problems of the lymphatic system usually cause swelling and other symptoms.

BTL's Lymphastim works on the pneumatic presso-therapy principle. Special applicator sleeves with multiple overlapping chambers provide a gentle massage that encourages the natural circulation of the lymph through the body. The effect is similar to the manual lymph drainage given by a massage therapist's hands.

## DIFFERENT FIELDS OF APPLICATION

- ◇ Rehabilitation and sports medicine
- ◇ Lymphology and rheumatology clinics
- ◇ Spas, wellness hotels and aesthetic/beauty centers

## LYMPH DRAINAGE ADVANTAGES

- ◇ Well-established & clinically-proven therapeutic method
- ◇ Immediate, visible results
- ◇ The treatment is well-received by patients
- ◇ Body is not burdened with unnecessary medication

## MEDICAL EFFECTS

Lymphatic drainage unblocks lymphatic ducts and vessels. This causes the natural flow of the lymph from the peripheral parts of the body, thus it reduces swelling and pain caused by the excessive amount of the lymph and metabolites in the limbs.

Lymphastim application also helps by healing and regenerating muscle tissue, increasing skin tone and restoring the original elasticity of the skin.


## MOST COMMON INDICATIONS

### Medical


Chronic oedemas  
Post-surgical obstruction of lymph vessels  
Varices prevention and treatment  
Venous insufficiency and ischemic diseases

### Wellness and aesthetics

Cellulite  
Postliposuction care  
Regeneration and detoxification  
Heavy-legs syndrome

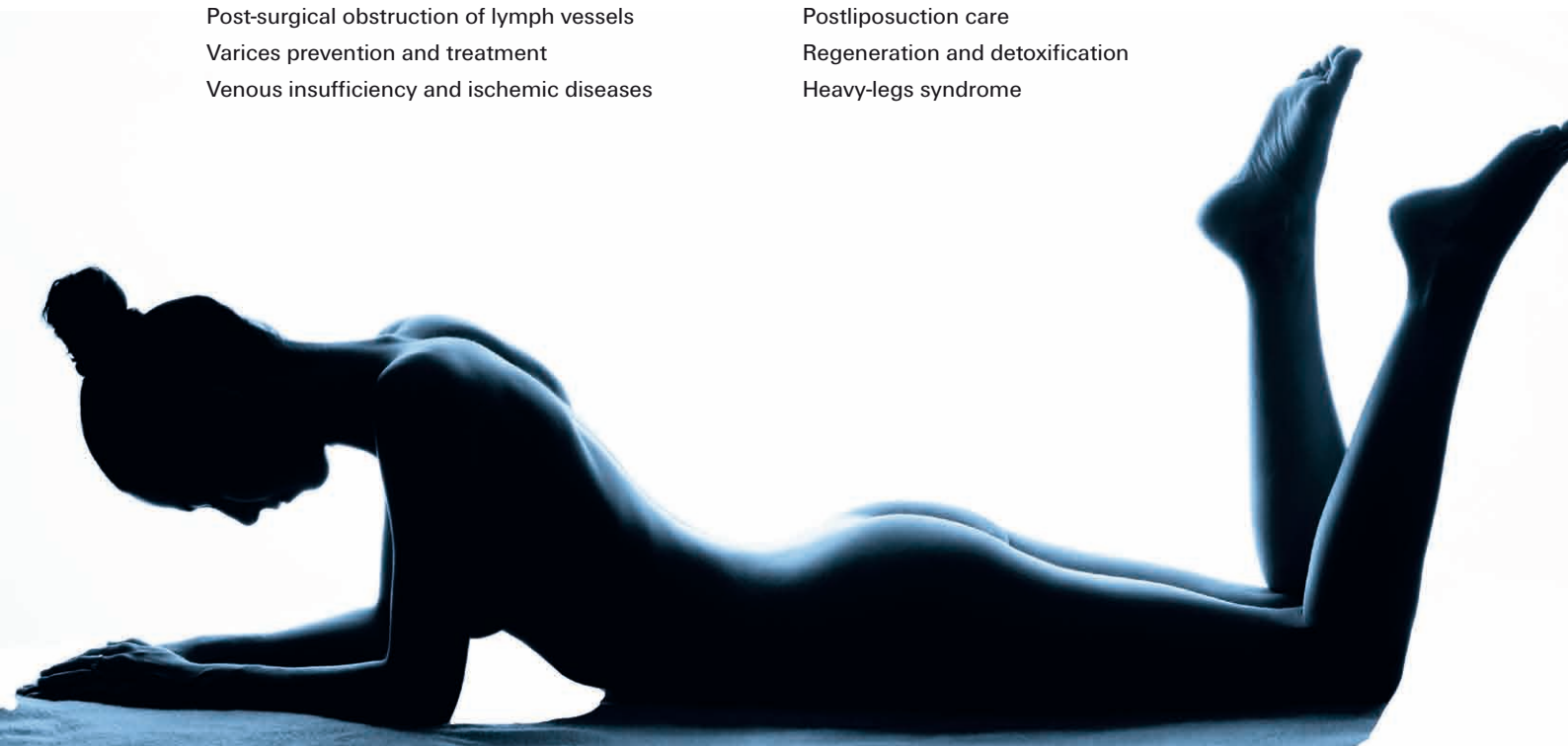

## HIGH QUALITY APPLICATORS

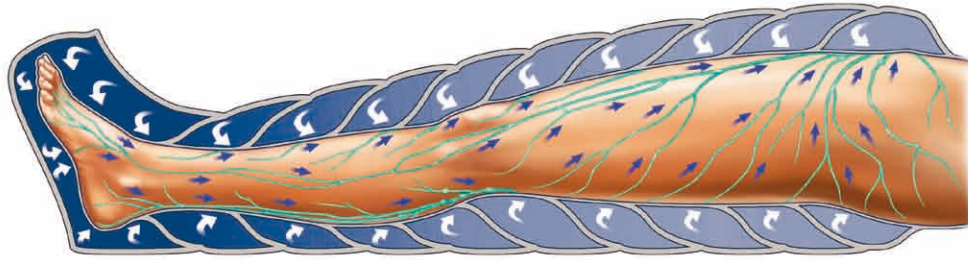

- ◇ **Fully overlapping sections**, without gaps, for a gentler and more effective movement of lymphatic fluids
- ◇ Extremely **abrasion resistant** polyamide material provides a unique combination of light-weight and durability
- ◇ **Antibacterial and dirt-repellent** material specially designed for use in a medical environment
- ◇ Ergonomically designed for maximum comfort

## EASY TO USE

- ◇ Connectors allow easy applicator swapping. Just plug in and treat.



- ◇ System has automatic identification of the applicators; the unit automatically recognizes which applicator is connected to it and then chooses the right compression flow.

## EASY TO HANDLE

- ◇ Full-length, high quality zipper for quick and easy dressing / undressing
- ◇ Automatic sleeve emptying after treatment for comfortable handling



## Available applicators:

- ◇ Trousers 22 chambers / Trousers 11 chambers



- ◇ Leg 10 chambers / Leg 5 chambers




- ◇ Arm 8 chambers / Arm 6 chambers



## BTL-6000 LYMPHASTIM TOPLINE

Unrivalled professional 12-channel presso-therapy device

- ◇ User-friendly interface with **colour touch-screen**

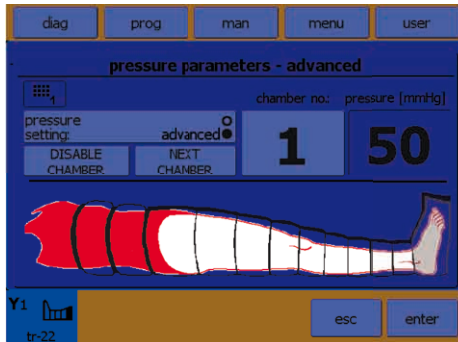

- ◇ **26 pre-set therapy protocols** and sequences of programmes. Protocols include unique pre-treatment programmes for opening of lymphatic nodes based on the principle of manual lymph drainage.
- ◇ 15 programs for various compression flow
- ◇ **Built-in encyclopaedia** with treatment recommendations
- ◇ Memory for up to 100 user-defined programmes
- ◇ **Advanced settings** for customized treatment include:
  - Individually adjustable pressure for each chamber
  - Working chamber activity shown on display
  - Option to disable one or more chambers where the patient cannot tolerate the pressure
  - Advanced settings of pressure gradient

## BTL-6000 LYMPHASTIM EASY

12-channel / 6-channel economy models with LED interface



- ◇ **Eight (8) pre-set programmes** that can be modified by the user
- ◇ Basic pressure gradient and programme speed setting
- ◇ Pressure can be changed during the course of therapy
- ◇ Working chamber activity shown on bar graph



**MODERN PRESSO-THERAPY SYSTEM  
WITH OUTSTANDING PARAMETERS**

## Technical Parameters



**BTL-6000**  
Lymphastim 12 TOPLINE



**BTL-6000**  
Lymphastim 12 EASY



**BTL-6000**  
Lymphastim 6 EASY

Ordering Number	P6000.003v100	P6000.002v100	P6000.001v100
User-interface	5.7" colour touch-screen	LED numbers	LED numbers
BTL Encyclopaedia	26 clinical protocols (programmes in sequences)		
User Programmes	100 user-defined programmes		
Pre-set Programmes	15	8	8
Number of Channels	12	12	6
Pressure Setting	20 -160 mmHg	20 -160 mmHg	20 -160 mmHg
Gradient	0 -100% continuous setting	0, 10, 20, 30%	0, 10, 20, 30%
Standard Accessories	Basic tube 12 TOPLINE, 1x power cord & touch-screen pen pointer	Basic tube 12 EASY & 1x power cord	Basic tube 6 EASY & 1x power cord
Optional Accessories	Trousers 22 chambers 2x leg + belt 2 chambers  Leg 10 chambers  Arm 8 chambers	Trousers 22 chambers 2x leg + belt 2 chambers  Leg 10 chambers  Arm 8 chambers	Trousers 11 chambers 2x leg + belt 1 chamber  Leg 5 chambers  Arm 6 chambers
	Interface for 2 arms or 2 legs: For using 2 arm or 2 leg applicators together		
	Extension strips (10 cm): For enlarging the circumference of the applicator		
	Trolley		



sales@btlnet.com  
www.btlnet.com

All rights reserved. Although every care has been taken to provide accurate and up-to-date information, no responsibility can be accepted for any mistakes or misprints contained in this catalogue. The products and their specifications, as shown in this catalogue, may change without notice.